

# Patient Mobility Products

*Godrej* interio

## Launch Circular for

- **EMERGENCY RECOVERY TROLLEY**
- **TRAUMACARE TROLLEY**
- **FOLDABLE WHEELCHAIR**
- **FOLDABLE WHEELCHAIR WITH CUSHION**





# Patient Mobility Products

When Patients arrive at the Hospital for emergency treatment , they are frequently brought by Ambulance, especially in cases of trauma / accidents / cardiac arrest / stroke / any condition which renders the patient unconscious / immobile.

In Emergency Cases, safe and swift mobilization of the patient is critical, from the Ambulance to the Operating Room and further from the Operating Room to the ICU / Ward .

For Paramedics, every minute counts, and the patient needs to be maintained in a stable condition even during mobilization. i.e. Blood, Glucose, and Oxygen supply need to be maintained steadily.

Further, during various stages of treatment and recovery, Patient Mobilization within the Hospital premises is done through stretchers and wheelchairs.

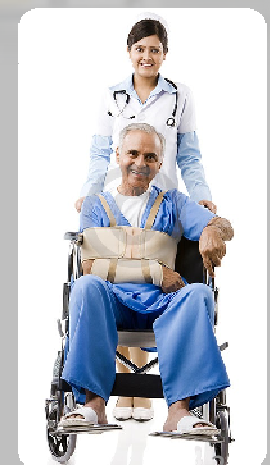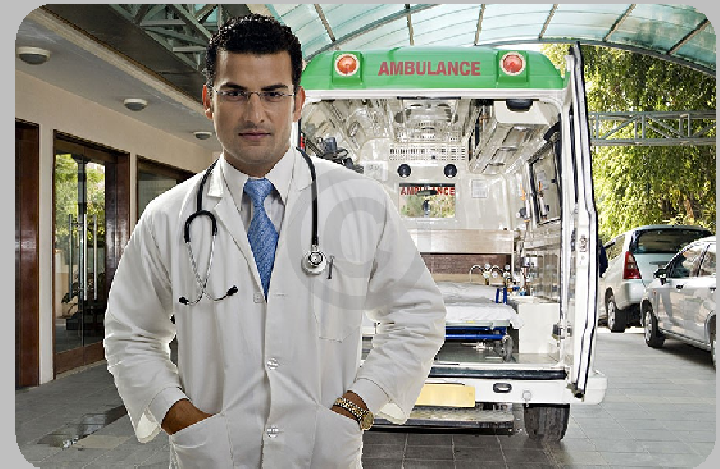

# Patient Mobility Products

Therefore the Traumacare and Emergency Trolleys used for this purpose are usually fitted with accessories like IV stands and Oxygen Cylinder Holders.

Further, during various stages of treatment and recovery, Patient Mobilization within the Hospital premises is done through stretchers and wheelchairs.

Accordingly, Wheelchairs may be used for Short trips to and from OPD or Ward, and later for Long Term use i.e. for Paralytic patients / Patients undergoing a long term recovery eg. Leg fractures



# Patient Mobility Products

## Some essential features

- **Patient Safety** : Considering that Patients are going to be moved sometimes at high speed, the product should shelter adequately from the sides to prevent patients from slipping / falling over / limbs getting caught. Safety bumpers are also preferable, to absorb shock in case of impact / collision.
- **Provision for multiple accessories**: The product should have provision for mounting of life-saving equipment like Oxygen Cylinders, IV Rods, etc.
- **Ease of mobility**: The product should have sturdy, well-aligned castors with high ground clearance for ease of mobility. If possible, propelling handles and handbrakes should be provided.

# Patient Mobility Products



## EMERGENCY & RECOVERY TROLLEY

- Hydraulic column type trolley with complete X-ray permeable surface.
- Overall Size: (L) 2130mm x (W) 800mm x (H) 560 to 910 mm
- Backrest and Kneerest adjustable on Gas Spring,
- Foot rest on Ratchet,
- Height Adjustment on Hydraulic Foot Pedal.
- High Coverage Side rails for enhanced safety.
- 5<sup>th</sup> Castor along with 200mm castors with directional lock provide ease of mobility.
- Push handle provided at front and rear side for the ease of trolley transfer with the better grip.
- ABS chassis ABS chassis at base provided with oxygen cylinder holder and store compartment.
- Max Patient Load: 135 kg

# Patient Mobility Products



## **TRAUMACARE TROLLEY**

- SS frame trolley with mechanical crank operations for Height Adjustment, Backrest and Knee rest.
- Overall Size: (L) 1900mm x (W) 750mm x (H) 520to 750 mm
- Backrest and Kneerest adjustable on crank
- Height Adjustment on Crank.
- High Coverage Side rails for enhanced safety.
- 150 mm dia castors with single castor with directional lock provide ease of mobility.
- Complete with Mattress, IV Pole holder, oxygen cylinder holder and storage shelf
- Max Patient Load: 135 kg

# Patient Mobility Products



## **FOLDABLE WHEELCHAIR WITH CUSHION**

- Wheelchair with 24inch rear mag wheels with alloy in the rim
- Front wheels are made of PA polymer & outer part is of solid rubber
- Arm rest are made of PU mould & base is made of ABS
- Seat is made of leatherite cushion seat which is rigid and durable, yet provided with a hinge in middle for folding.
- Hand brakes are provided
- Handles are made of molded rubber for better grip

# Patient Mobility Products



## **FOLDABLE WHEELCHAIR**

- Wheelchair with 24inch rear mag wheels with alloy in the rim
- Front wheels are made of PA polymer & outer part is of solid rubber
- Arm rest are made of PU mould & base is made of ABS
- Seat is made of leatherite material seat which is flexible and foldable.
- Hand brakes are provided
- Handles are made of molded rubber for better grip



# Patient Mobility Products

SR #	MODEL	CPL [Rs]	WPL[Rs]	ITEM GROUP	ITEM CODE
1	TRAUMA CARE TROLLEY	116350	95407	DTHOTS	STHC5X4100000000
2	EMERGENCY RECOVERY TROLLEY	288000	236160	DTHOTS	STHC5X5200000000
3	FOLDABLE WHEEL CHAIR WITH CUSHION	12500	10250	DTHWHC	STHC4XA221000001
4	FOLDABLE WHEEL CHAIR	12500	10250	DTHWHC	STHC4XA211000001

**NOTE:**

*The given CPL & WPL values are NOT inclusive of Taxes. Please include Tax/other expenses over & above the given values as applicable.*

*ED is Nil.*

*Delivery Terms are 4-6 weeks from Booking of Order.*

# Patient Mobility Products

- For Product Related Technical Queries / Sales Support, please contact Swati / Samrat at 9920352879/ 9833729210 & at [swatiw@godrej.com](mailto:swatiw@godrej.com)/ [samratg@godrej.com](mailto:samratg@godrej.com)
- For Order Booking, please contact Anshuman / Ruchi at 9833558260 / 9969820899 & at [ams@godrej.com](mailto:ams@godrej.com)/ [ruchi@godrej.com](mailto:ruchi@godrej.com)
- For Dispatches, please contact Joseph at 9819940733 & at [jjd@godrej.com](mailto:jjd@godrej.com)
- For Product Specifications in pdf, please check Intranet after a week.
- For APOs, please log entry in the COTS with relevant details.

**Good Luck & Happy Selling!**